


SPECIFICATIONS:

HOSE ASSEMBLIES

The Hose Assemblies shall be an Amico **Alert-1** series.

The Hose Assemblies are manufactured with FDA approved medical grade thermoplastic conductive rubber tubing. The tubing is reinforced with nylon braiding. The maximum working pressure is 200 psi @ 70° F. The tubing is color coded to conform to CGA C-9 standard. The fittings are permanently crimped on each end with ferrules.

The Hose Assemblies are normally made with a DISS Female Handtight nut on one end (to be connected to an outlet). The other end is a Male DISS or Female Quick Connect coupler with Check Unit. The length is measured in feet from end to end (including fittings).

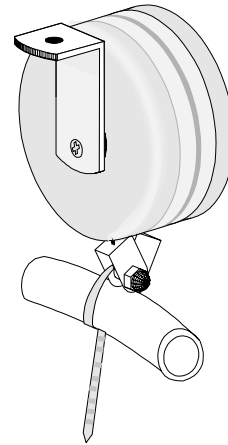

SPECIFICATIONS:

HOSE RETRACTOR

The Hose Retractor shall be an Amico **Alert-1** series.

The Hose Retractor can be attached directly to a ceiling outlet or directly to the ceiling. The retractor wire is manufactured out of 48" stainless steel heavy duty cable.

MODEL NUMBER:

SHC-05C-DFOF-A

Length in feet:



Color code:

- U = US
- E = ISO

End Connections:

- Chemetron = C
- DISS = D
- Ohmeda = O
- Puritan = P
- NPT = N
- Male = M
- Female = F

Gas:

- Oxygen = O
- Medical Air = A
- MedVac = V
- Nitrous Oxide = 2
- Nitrogen = N
- Carbon Dioxide = C
- WAGD/AGSS = E

For example an Oxygen 5' (60") Hose Assembly with Female DISS by Ohmeda Female coupler with US Colours:

SHC-05U-DFOF-O

MODEL NUMBER:

S-HOSE-RET-HD

PROJECT:

D I S T R I B U T E D B Y :